

ULTRA TRANSLUCENT

UZir is made from a unique formula of Ultra Translucent Zirconia.

CAD/CAM
milled
restorations

800 MPa
flexural
strength

scdlab.ie/UZir


Its enhanced translucency and strength allow you greater flexibility to create a natural and highly aesthetic result for your patient.

UZir is an ideal choice for anterior and posterior crowns as well as bridges of up to 3 units.


Fabricated using CAD/CAM technology, you can be assured of a precise fit and predictable result.

This innovative product is highly versatile and provides a wide range of indicators to meet your patient's needs.

Crown & Bridge	Single crowns and 3-unit bridges
Inlay/Onlay	Yes
Flexural strength	800 MPa
Stump Shade	Must be provided as it may affect the final shade. UZir is not recommended for cases where metal cores exist.
Cementation	Self-adhesive or adhesive resin cements are recommended. Shade of the cement will contribute to the final shade. Please contact our technical team for further information.
Adjustments	Fine diamond with water cooling or diamond impregnated rubber is recommended. Areas must be polished to a high shine with diamond polishing paste afterwards.

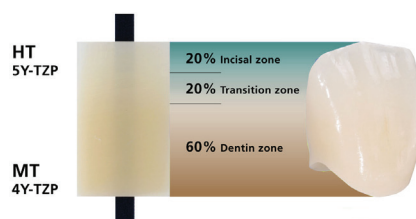

Translucency

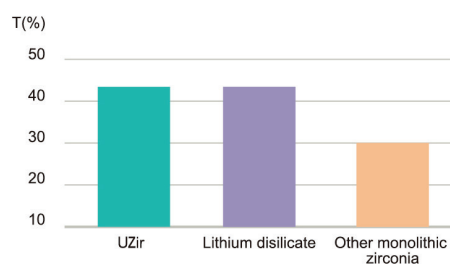

Flexural Strength

Three-point bending test

